

# Manufacturer/Supplier:

WUXI DONGLIN SCI&TECH DEVELOPMENT CO., LTD

Address: A1-203 Mingpin CityII, No.8 Xihudong Road. Liangxi District

Wuxi Jiangsu Province, China

TELE: 86-510-82732223 FAX: 86-510-82720101-8014

Web: www.dldevelop.com

Email: <u>info@dldevelop.com.cn</u> <u>service@dldevelop.com.cn</u>

**Product name:** Rabbit anti-Anamorsin Antibody

Catalog: DL96760A

**Synonyms:** Anamorsin; Cytokine-induced apoptosis inhibitor 1;

Fe-S cluster assembly protein DRE2 homolog

**Immunogen:** Recombinant full length protein of human

Anamorsin

Form: Liquid Concentration: 1mg/mL

Size: 100 ul/50 ul Host: Rabbit

**Reactivity:** Human, Mouse Application: WB, IHC, IF/IC

Clonality: Polyclonal Dilution: WB (1/500 - 1/2000),

IHC (1/50 - 1/200), IF/IC (1/50 -



#### 1/200)

**Entrez Gene:** 57019/109006

O6FI81/O8WTY4 SwissProt:

The antibody was purified by immunogen affinity **Purification:** 

chromatography.


**Buffer:** Liquid in 0.42% Potassium phosphate, 0.87%

Sodium chloride, pH 7.3, 30% glycerol, and 0.01%

sodium azide.

### **WB** description:

Western blot analysis of Anamorsin expression in HepG2 (A), Hela (B), MCF7 (C), mouse testis (D) whole cell lysates.

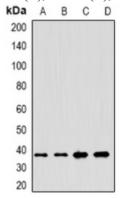


# **IHC description:**



Immunohistochemical analysis of Anamorsin staining in human liver cancer formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

#### **IF/ICC description:**

Immunofluorescent analysis of Anamorsin staining in HepG2 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary anti

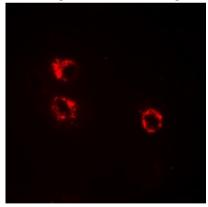

### **Storage:**

Store at -20°C. Avoid repeated freeze / thaw cycles.