

Manufacturer/Supplier:

WUXI DONGLIN SCI&TECH DEVELOPMENT CO., LTD

Address: A1-203 Mingpin CityII, No.8 Xihudong Road. Liangxi District

Wuxi Jiangsu Province, China

TELE: 86-510-82732223 FAX: 86-510-82720101-8014

Web: www.dldevelop.com

Email: <u>info@dldevelop.com.cn</u> <u>service@dldevelop.com.cn</u>

Product name: Rabbit anti-DDX58 Antibody

Catalog: DL96792A

Synonyms: Probable ATP-dependent RNA helicase DDX58;

DEAD box protein 58; RIG-I-like receptor 1; RLR-1; Retinoic acid-inducible gene 1 protein; RIG-1; Retinoic acid-inducible gene I protein;

RIG-I

Immunogen: KLH-conjugated synthetic peptide of human

DDX58

Form: Liquid Concentration: 1mg/mL

Size: 100 ul/50 ul Host: Rabbit

Reactivity: Human Application: WB, IHC, IF/IC



Clonality: Polyclonal **Dilution:** WB (1/500 - 1/2000),

IHC (1/50 - 1/200), IF/IC (1/20 -

1/100)

Entrez Gene: 23586

O95786 SwissProt:

Purification: The antibody was purified by immunogen affinity

chromatography.

Buffer: Liquid in 0.42% Potassium phosphate, 0.87%


Sodium chloride, pH 7.3, 30% glycerol, and 0.01%

sodium azide.

WB description:

Western blot analysis of DDX58 expression in Jurkat (A), HepG2 (B), MCF7 (C) whole cell lysates.



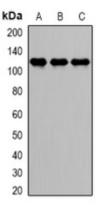


IHC description:

Immunohistochemical analysis of DDX58 staining in human prostate formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

IF/ICC description:

Immunofluorescent analysis of DDX58 staining in HepG2 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody



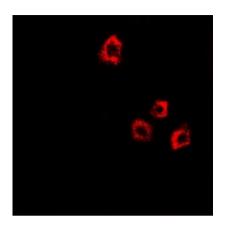

Storage: Store at -20°C. Avoid repeated freeze / thaw cycles.